



Endless  
**LifeStyle**  
Possibilities



Maximize  
Your  
**Oxygen**  
Mobility...



...Like Never Before

LIFESTYLE™ PORTABLE OXYGEN CONCENTRATOR  
FOR ACTIVE, WORKING, AND  
TRAVELING OXYGEN PATIENTS



# LifeStyle™ by AirSep®

Oxygen patients and their caregivers all over the world are now traveling lighter, further from home, and free of the restrictions formerly imposed by oxygen refills and deliveries.



Whether it's a trip to the market, an unplanned overnight or extended stay with relatives, or just out for coffee with a dear friend, no other oxygen system can meet your needs better or more conveniently than LifeStyle.

## Portable Oxygen Solution

LifeStyle by AirSep is the consummate portable oxygen solution — a less than 10 lb. (4.4 kg) concentrator that actually makes oxygen everywhere you go. AirSep is a leading manufacturer of high-quality respiratory equipment for home, hospital, and clinical use.

If you are an active oxygen patient, you may have already experienced the limitations of needing to bring all of your contents with you, and still occasionally felt the anxiety of running out of oxygen due to circumstances beyond your control.



In contrast, you will feel liberated by LifeStyle's unique oxygen producing capability, which is truly irreplaceable. This unique system enables you to return to work, join the family reunion, fly via commercial aircraft, or travel on an exotic cruise to meet your heart's desire.

## No Cylinders, No Refills – Ever Again!

Imagine – no need for cylinders or liquid system refills! That means no coordinating deliveries at home or away and no patiently awaiting a visit from the equipment provider to ensure your contents are filled so that you're ready to go — with LifeStyle you already are!

## Multiple Power Options

Multiple power options mean you can be outdoors or just about anywhere indefinitely with simple battery recharges or replacements. The unit operates continuously and recharges the




battery on normal household AC electric or from the 12-volt DC outlet available in most automobiles or motorized boats.

Thousands of patients around the globe now use the one-of-a-kind LifeStyle POC (Portable Oxygen Concentrator)

to meet their daily supplemental oxygen needs. And LifeStyle is durable enough to use 24 hours a day.

Working patients, traveling patients, and those for whom the active lifestyle is vitally important rely on the safety, convenience, and versatility of the most modern and technologically advanced method of oxygen delivery — via POCs.

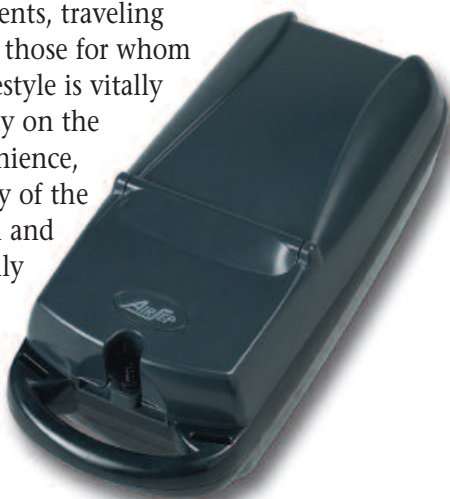



**Rechargeable Battery**  
BT007-1

**AC Power Supply**  
PW005-1

**Power Cord**  
CD023-2 (120 Volt hospital-grade)  
CD017-2 (220 Volt)  
CD025-1 (220 Volt Australia)

**DC Auto Adapter**  
PW015-1

**LifeStyle Retractable Cart**  
MI210-1

LifeStyle Portable Oxygen Concentrator		PowerPack
AS081-1	Complete with 120V power cord	BT011-1
AS081-3	Complete with 220V power cord	BT011-3
AS081-4	Complete with 220V power cord – Australia	BT011-4

## SPECIFICATIONS

**OXYGEN CONCENTRATION\*:** 1-5 pulse settings; equivalent to a continuous flow of 90% oxygen at  $\pm 3$

**DIMENSIONS:** 5.5 in. H x 7.25 in. W x 16.31 in. L  
(13.9 cm H x 18.4 cm W x 41.4 cm L)

**WEIGHT:** 9.75 lb (4.4 kg)

**POWER:** LifeStyle unit — 13.5 VDC, 5.5 amps  
Power supply — 100-250 VAC,  
50/60 Hz, 1.0 amps maximum

Battery pack — 12 VDC, nominal

**BATTERY DURATION:** 50 minutes nominal for fully-charged battery

**BATTERY QUICK-CHARGE:** 2 to 2-1/2 hours

**WARRANTY:** 2-year parts on unit;  
90 days on rechargeable battery

\*Based on atmospheric pressure of 14.7 psia (101 kPa) at 70° F (21° C)

Using efficient pulse-flow technology, LifeStyle delivers oxygen only on inspiration, thereby eliminating the wastefulness incurred by continuous flow systems. This advanced technology makes possible this lightweight, compact, portable concentrator to accompany you anywhere opportunity takes you.

## Lightweight and Ready for Adventure

The first truly portable oxygen concentrator of its kind — LifeStyle — has encouraged patients, loved ones, and other caregivers to do what they have not done in years, or to simply try something new, all with maximum mobility.



## Easy to Operate

This easy-to-use unit allows you to be up and running in no time. Just flip open LifeStyle's cover, attach the end of your oxygen cannula to LifeStyle's oxygen outlet, turn the unit on to your prescribed setting, and you're all ready to go. What could be easier?

## Alarm System

The unit's intermittent green light display shows that you are receiving a pulse flow on every breath. The unit also illuminates with an intermittent yellow indicator and then audibly alarms to alert you when the battery is near depletion. The red indicator, which can be accompanied by an audible alarm as well, signals that the cannula may be inadvertently disconnected, there is a kink in the tubing, or no breathing is otherwise detected by the unit at that time — all valuable safety features.

## A Custom-Made Cart for Added Convenience

A portable, retractable cart gives you even more convenience and flexibility. There's no need to

carry the unit when out for a short stroll or when taking an extended walk.

The cart also features a custom-rack system, which holds PowerPack™ a multi-battery/charger.



## The Perfect Accompaniment for Airport Logistics

LifeStyle eliminates the difficulties of air travel. Think about never again needing to plan airport deliveries of oxygen. Schedule and make your flight connections with LifeStyle instead, and there's no need to panic about unexpected delays or layovers ever again.

Note: As the first POC, LifeStyle was also one of the first to receive the US Federal Aviation Administration's (FAA) acceptance for onboard, in-flight use by oxygen passengers on commercial aircraft for effectivity August 11, 2005 via SFAR 106.

For a list of airlines that accept passengers with LifeStyle POCs, reference [www.airsep.com](http://www.airsep.com). Please check directly with the individual airlines with which you would like to travel for up-to-date information on their specific POC policies.

To upgrade or to check on the status of your POC for air travel, contact your oxygen equipment provider.



## For Some, a Luxury; For Others, a Modern Medical Convenience



At this point in your life, convenience and greater mobility may further enrich the quantity and quality of the time you and your loved ones are able to spend away from home. Whether it's simply a trip to the beauty salon, an extended holiday or

vacation, such as a special anniversary cruise, or the travel back and forth to a second home, LifeStyle can fit right in with all your journeys.

## For the Working Patient


For those who wish to maintain a regular work schedule away from home despite the need for



supplemental oxygen, LifeStyle may again serve as the perfect choice. Drive back and forth to work with LifeStyle connected to your automobile's standard 12-volt DC connection, use the unit's internal battery when leaving the car, or going between

offices and other buildings, and then connect into any standard AC electrical outlet nearby.

The AC or DC power automatically recharges the internal battery while the unit remains operational, or you can replace a depleted battery with a fully charged one and be off on your way again.

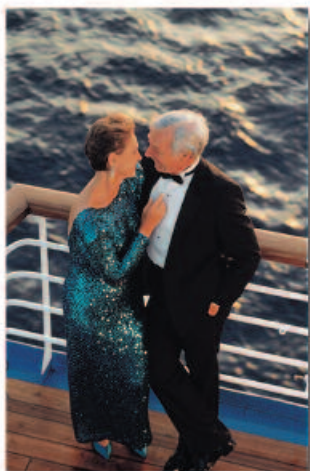

The beauty of LifeStyle is that wherever, whenever you go – you're never away from a power source and its endless supply. LifeStyle continues to make your oxygen contents, delivering just the amount of supplemental oxygen your prescription requires, as you need it.

## Low Power Consumption

At just 42 watts, LifeStyle's low power consumption is an added plus, for home or clinical use, as it draws less power than an average light bulb, and merely a fraction of the power consumed by the average full-size concentrator.




## Use in Clinical Settings

For those in nursing homes, assisted living, or other institutional settings, LifeStyle is, again, the perfect accompaniment. Stroll the grounds, move through the dining halls or recreational centers, alone or with family and friends, and



discover the true freedom of LifeStyle. Even a trip away is now made easier. Without hesitation, you know that you and your oxygen are always ready to go -- off to greet the adventures of a new day.

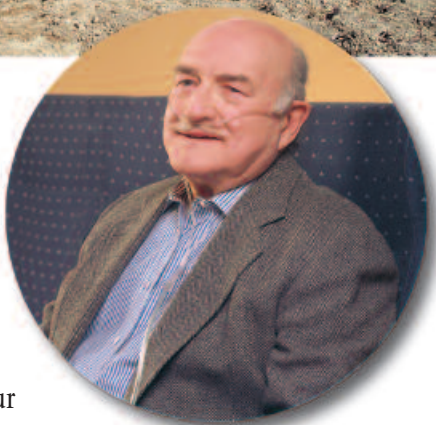

## **Works Well for Remote Areas**

If you're one of the hundreds of thousands of supplemental oxygen users residing in remote areas, or are not able to receive deliveries from an equipment provider frequently, this LifeStyle choice becomes extremely valuable. Don't feel rationed or go without oxygen if you have been prescribed for it. Proper oxygen levels are a necessity for maintaining your current health and well being and may extend your life expectancy, particularly in conjunction with daily ambulation. LifeStyle ensures that you are free to ambulate for as often and as long as you feel you can.

## At Rest or At Play

You, your physician and equipment provider will determine if LifeStyle is suitable for your specific oxygen condition. And, as with any conserving device,

AirSep recommends the use of pulse oximetry to determine the appropriate setting(s) that ensure the proper oxygen saturation levels for individuals at rest, during exercise or exertion, or even during sleep, in accordance with health care professionals' guidelines.



LifeStyle's "cannula disconnect" alarm will provide a continuous alarm should the cannula become disconnected inadvertently (no inspiration detected), or if the unit is not being used while it is running.

For certain patients (such as mouth-breathers), it may even be desirable to wear a chin strap.

## How Do I Pay for LifeStyle?

LifeStyle is exclusively marketed through authorized oxygen equipment providers; therefore, your local provider may offer a variety of programs for the purchase, or short- or long-term rental of this product. Depending on your insurance policy and coverage, there may also be full or partial reimbursement for LifeStyle. Discuss with your physician, equipment provider, and insurance company today the possibility of using the LifeStyle Portable Oxygen Concentrator to ensure a whole new lifestyle tomorrow.

*P*icture the future with AirSep's LifeStyle and you at the center of it. Start today to plan your future. Dreams do become reality when you reach for new ways to access your lifestyle choices.



**Don't delay!**

*C*ontact your equipment provider to find out how you can purchase or rent LifeStyle for you or a loved one's oxygen needs.



AirSep Corporation  
401 Creekside Drive  
Buffalo, NY 14228-2085 USA